Water Purifier - Auto Pure Type 1 Water

For Genetics, DNA Sequencing, Polymerase Chain Reaction (PCR)

Benchtop WG301UVUFB: 120V / WG311UVUFB: 220V Remote Dispense WG301UVUFR: 120V / WG311UVUFR: 220V

Water quality Type 1, 18.2 M Ω

Feed Requirement < 20 uS

Flow

0.67 gpm 2.5 lpm

Water is purified using a staged purification process which includes high-purity ion exchange resins to remove dissolved minerals, ultraviolet light for bacteria sterilization and TOC reduction, and ultrafiltration for pyrogen removal and nuclease free water. As water exits the system, a final filter removes particulates and bacteria to reach Type I water specifications.



■ WG Series Configuration

- Conditioning cartridge
- Polishing cartridge
- UV oxidation
- Ultrafiltration
- 0.2 micron filter


Features

- 2.5 LPM* of 18 MΩ Type 1 reagent water on demand
- Intuitive, user-friendly touch screen interface
- Real time quality monitoring with leak detection
- Compact design with remote dispense option
- Quick change, no tools cartridge design
- Large capacity, low cost consumable cartridges
- Premium quality ion exchange resin
- Low ownership cost

Specifications

Model	WG301UVUFB/WG311UVUFB	WG301UVUFR/WG311UVUFR
Туре	Benchtop	Remote dispense
Resistivity*	18.2 MΩ-cm	
Bacteria*	< 1 cfu/ml	
Endotoxin**	<0.005 EU/m	
Particulates*	< 0.05 µm filtration	
TOC*	< 5 ppb	
RNase*	< 0.01 ng/ml	
Dnase*	< 4 pg/µl	
Temperature	100°F / 30°C	
Pressure	90 PSIG Max. / 20 PSIG Min.	
Flow rate (typical)	0.67 Gallons Per Minute / 2.5 Liters Per Minute , < 2.0 LPM with Ultrafilter & Endotoxin filter installed	
Inlet connection	3/8" O.D. Tubing / 9.53 mm O.D. Tubing	
Dimensions (H x W x D)	21.6 in. x 14.7 in. x 14.4 in. (54.8 cm x 37 cm x 36.5 cm)	
Weight	42 lbs. dry / 47 lbs. operating (22.7 kg dry / 25 kg operating)	
Power requirements	120V, 60 Hz @ 1.0 amp / 220V, 50 Hz @ 0.5 amp	

^{*} Performance based upon RO / DI or Distillation Feed (< 20 uS/cm, < 50 ppb of TOC, 0.2 micron filtration) ** Includes Endotoxin Removal Capsule Filter

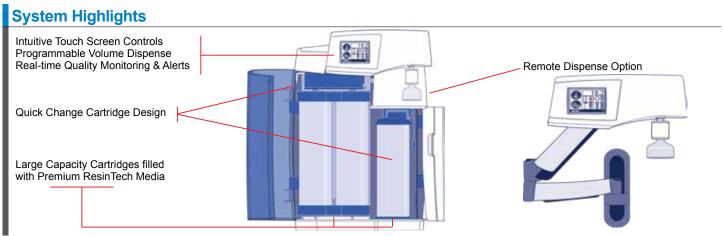

^{*} Reported flow rate is typical but can vary depending on supply pressure and system options

Flow Diagram Direction of Flow TOC Monitor Solenoid Valve (Optional) Ultrafilter Pump Check Valve 0.2 micron Filtration Intermediate Resistivity Sensor Final Resistivity Sensor Inlet Pressure Regulating Valve Endotoxin Filter Final Filter Conditioning 185/254 nm Polishing Optional Flow Sensor Cartridge Cartridge UV Lamp

Additional System Components

Description		
Installed Options and Accessories*		
Direct Feed Port		
Recirculating Dispensing Gun		
Wall Mount Bracket		
Cartridges and Filters**		
Dual Conditioning Cartridge (RO Feed)		
Dual Conditioning Cartridge (SDI / TAP Feed)		
Polishing Ultrapure Cartridge		
Polishing Ultrapure Low Organics Cartridge		
0.2 Micron Capsule Final Filter - 1/4" MNPT x 1/4" hose barb		
0.2 Micron Capsule Final Endotoxin Filter - 1/4" MNPT x 1/4" hose barb		
Consumables		
Ultrafilter - 3/8" Quick Disconnect Fitting		
UV Bulb 254/185nm TOC Destruct		
Sanitization Kit (includes cartridges, test strips, and gloves)		

^{*} Must be pre-ordered for factory installation

■ Other WG Series Configuration (includes TOC Monitor)

Product code and Description	
WG301UVUFTB	
Benchtop - with UV, UF and TOC Monitor installed, 120V/60Hz	
WG311UVUFTB	
Benchtop - with UV, UF and TOC Monitor installed, 220V/50Hz	
WG301UVUFTR	
Remote Dispense Configuration - with UV, UF and TOC Monitor installed, 120V/60Hz	
WG311UVUFTR	
Remote Dispense Configuration - with UV, UF and TOC Monitor installed, 220V/50Hz	

Best Value in Replacement Consumables



^{**} Not included in a system and orderered separately